


Peptide Synthesis Report

MOFFITT
CANCER CENTER

PROTEOMICS

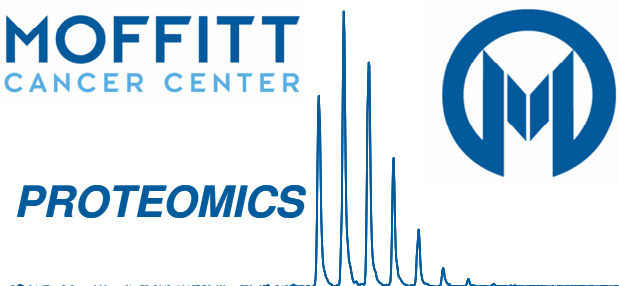

Project Title: Internal standard for Casein kinase I isoform alpha (KC1A_HUMAN)

Submitter: Yi Chen

PI: John Koomen

Design and Synthesis:


MPS#	Peptide Sequence	Residues	Modification	MWT Monoisotopic
p-0111	LFLIDFGLGK	153-162	A9 to G	1121.64852

Solid state peptide synthesis (Symphony, Protein Technologies, Tucson, AZ) is used to make standards at the 25 micromole scale using standard FMOC chemistry. Briefly, Wang resin (Novabiochem, Germany) serves as the support to initiate peptide formation. Activator solvents for the amino acid coupling consist of 150 mM N-Methylmorpholine in N-methylpyrrolidone (NMP) with 100 mM of O-benzotriazol-1-yl-N,N,N',N'-tetramethyluronium hexafluorophosphate (HBTU). Deprotection between couplings is achieved using 20% piperidine in NMP with 2% 1,8-diazabicyclo [5,4,0] undec-7-ene. Upon sequence completion, cleavage from the resin is performed using 95% TFA, 2.5% water, and 2.5% triisopropylsilane. Peptides are then precipitated in ice cold ethyl ether, washed twice, and resuspended in water prior to lyophilization.

Difficult, intermediate, easy synthesis coupling map:
LFLIDFGLAK to LFLIDFGLGK

Semi-preparative HPLC:


HPLC purification is performed on a semi-preparative system (U3000, Dionex, Sunnyvale, CA) using a C18 reverse phase column (TP 238, 250 x 10mm, 10-15 μ m particles, Grace Vydac, Deerfield, IL). Twenty minute gradients are run from 5% to 60% B solvent (A: 2% acetonitrile/0.1% formic acid; B: 100% acetonitrile/0.1% formic acid). Eluted peptides are measured at 214 nanometer wavelength and collected with automatic triggering (Foxy Jr, ISCO).

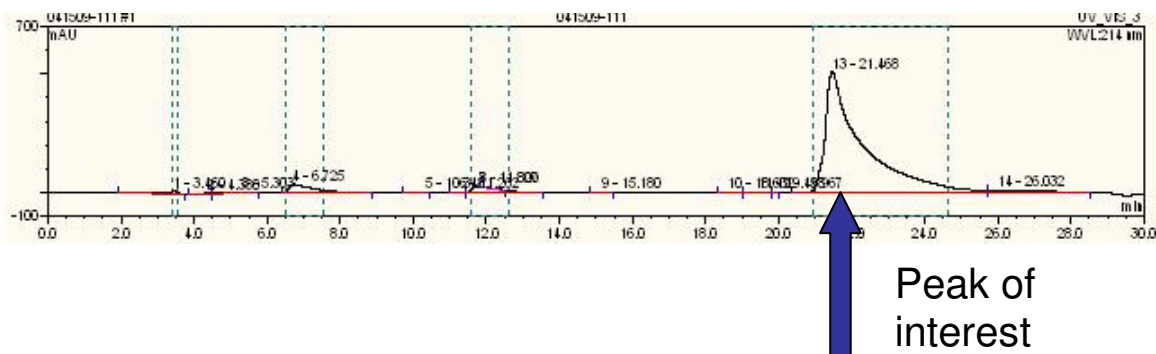


MALDI:

Purified peptides are analyzed with MALDI MS to verify purity and sequenced with MS/MS (4700, Applied Biosystems, Framingham, MA). Peptides in 2% ACN with 0.1% formic acid are mixed 1:1 with α -cyano-4-hydroxycinnamic acid (10 mg/ml) in 50% ACN and deposited in 1 microliter aliquots on the MALDI target.

MS Spectrum: ~100% purity

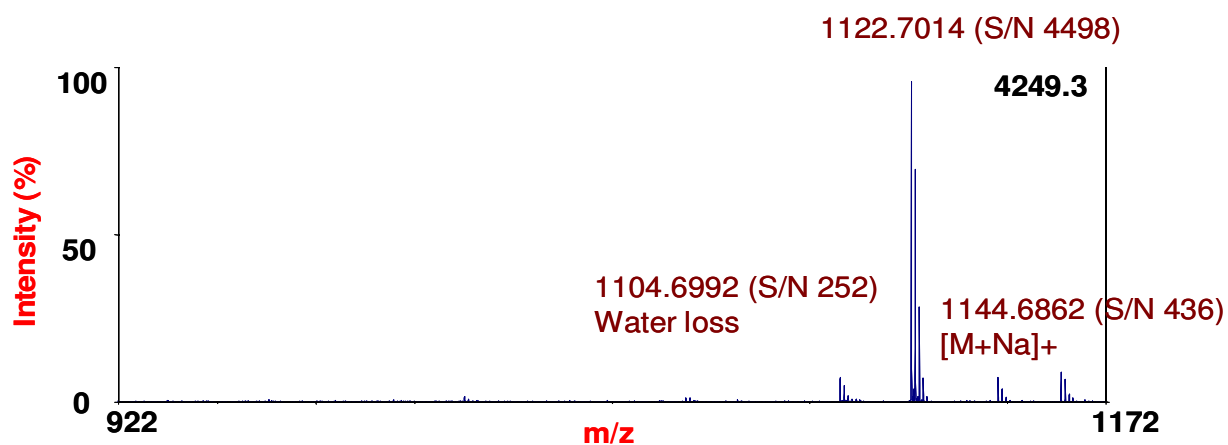

MS/MS spectrum:

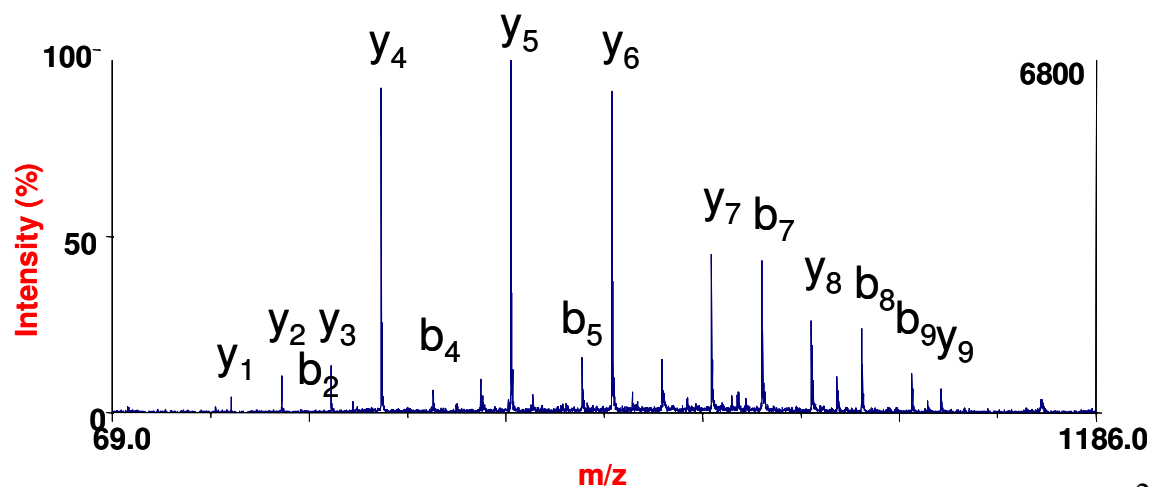

Fragment Ion table with Identified peaks:

Seq	#	B	Y	#
L	1	114.092	1122.656	10
F	2	261.160	1009.572	9
L	3	374.244	862.504	8
I	4	487.328	749.420	7
D	5	602.355	636.336	6
F	6	749.424	521.309	5
G	7	806.445	374.240	4
L	8	919.529	317.219	3
G	9	976.551	204.135	2
K	10	1104.646	147.113	1